




Jordco® e-Ruler®

Endodontic File Measuring Ruler
Available in Two Formats and
Two Colors

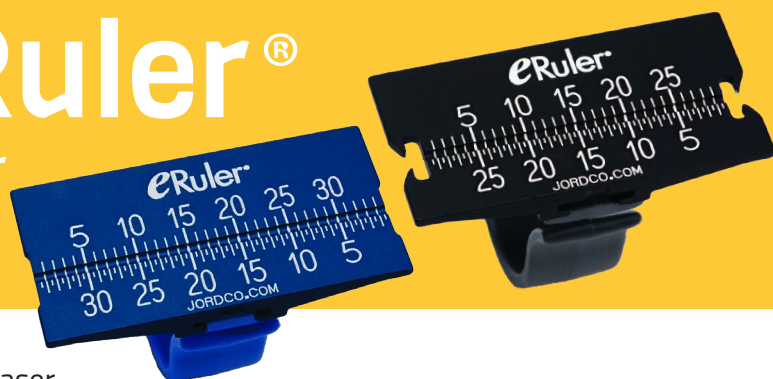

The laser etched e-Ruler instrument is worn on the clinician's finger. The ruler's markings are precisely laser etched providing enhanced contrast for accurate and easy calibration of endodontic hand and rotary files. The New REDUCED GLARE black color reduces head lamp glare while operating with surgical telescopes.

TECHNICAL SUMMARY

e-Ruler: The e-Ruler instrument is constructed of blue anodized aluminum with a permanent laser etched scale. It comes with three sizes of flexible thermoplastic rings that can be tried for fit prior to snapping the correct ring size into the e-Ruler instrument. Two versions of the e-Ruler instrument are available including a 30 mm ruler with stop locks and a 35 mm ruler with no stop locks.

EndoRing II Endodontic Instrument Family of Products

The EndoRing II organizer enables the clinician to directly place, store, measure and clean endodontic hand and rotary instruments within the operating field. It reduces the risk of cross-contamination and minimizes the passing of sharps between dentists and staff. Using the EndoRing II organizer helps make endodontic procedures safer and faster.



The EndoRing FileCaddy® organizer provides a unique bulk storage solution for all endodontic files. The clinician places multiple files securely into the special foam insert, closes the lid and steam autoclave processes the entire EndoRing FileCaddy organizer.



Discover our complete line of quality Jordco® products at www.jordco.com or contact your dental supply dealer.



Jordco, Incorporated
8531 NE Cornell Rd, STE 800, Hillsboro, Oregon 97124 USA
(800) 752-2812 toll free USA www.jordco.com
Visit our website for patent and trademark information.



PRODUCT SPECS

Name: e-Ruler Pack with or without Stop Locks
REF: ERUL-s: Blue with Stop Locks
ERULX-s: Blue without Stop Locks
ERULRG-s: Black Reduced Glare with Stop Locks
ERULXRG-s: Black Reduced Glare without Stop Locks
Weight: .5 oz

PACKAGE CONTENTS

- 1 Metal e-Ruler
- 3 Finger Rings (S, M, L)

FEATURES AND BENEFITS

Laser Etched Scale

The e-Ruler scale is laser etched onto a blue anodized platform with precise, distinct millimeter and numeric markings. Both e-Ruler models include ½ mm marking on the lower half of the scale for extra fine calibration.

Rubber Stop Locks

Machined rubber stop locks are available on the 30 mm e-Ruler model. Rubber stops on the files are inserted into the rubber stop lock and the file can then be precisely calibrated using a simple push-pull technique.

Dual Right and Left Handed Ruler

Both e-Ruler models have laser etched numeric markings for use by either right or left handed operators.

Finger Held Ring

Comfortably worn on the finger, the e-Ruler allows precise measurement of hand and rotary instruments in the working field eliminating the need for measuring blocks or similar devices.

Easy Calibration of Curved Instruments

To measure a curved file, insert the instrument rubber stop into the e-Ruler rubber stop lock, align the instrument tip into the recessed calibrating channel and apply a slight downward pressure to flatten the file into the channel. The measurement is set using a simple push-pull technique and the instrument will return to its curved shape when removed from the e-Ruler.